

IUCLID 5.2 Helpsystem Plugin Installation Guide

© European Chemicals Agency, 2010

Published 10 June 2010

1. Introduction

This document provides a description of the installation procedure of the IUCLID 5 translated helpsystem plugins. More than one helpsystem plugin can be installed on one instance of IUCLID 5, however it is recommended not to install more than one language in addition to English, due to its resource requirements (memory and startup time).

The installation of each translated helpsystem consists of copying two files into the appropriate folders. The directory names may be different on your system depending on the operating system and version, but the concept is the same.



Note

In this document we use the French version of the helpsystem as an example. Please replace the two letter language code in the file names to the code of your own language.

2. How to Download the IUCLID 5 Helpsystem Plugins

- Navigate with your web browser to the IUCLID website [<http://iuclid.echa.europa.eu/>] and log in to your previously created account.
- In the “Download” section select “Version 5.2.0”, then “Helpsystem plugin” and click “Download the installation package”.
- Select the “I acknowledge” check box and then click the “Next” button.
- In order to download the software you have to accept the license agreement. Read carefully the License Agreement and, if you agree, accept it selecting both the “I accept” and “I am entitled to sign this License agreement” check boxes. Click the “Next” button. Select “please click here” to start the download. Select “Save” to save the file to your computer.
- Save the file `iuclid5_plugin_helpsystem_2010-02-12_fr.zip` to your computer (replace the two letter language code).
- Extract the contents of the downloaded file into a temporary folder:
 - `plugins/client/eu.echa.iuclid.plugin.helpsystemFR.zip`
 - `/plugins/server/eu.echa.iuclid.plugin.resourcesFR.zip`



Important

Do not unzip the two files above.

3. Installation Instructions



Note

Pre-condition: IUCLID release 5.2 installed on your computer.

3.1. Standalone Version

- If IUCLID 5 is running, save your work and close the application.
- Locate the installation folder of IUCLID 5 on your hard disk.




Note

The installation folder is usually C:\Program files\IUCLID5\5.2.0\, however it may be different on your system.

- Copy the file eu.echa.iuclid.plugin.resourcesFR.zip to the subfolder \plugins\server of IUCLID 5.


For example:

C:\Program files\IUCLID5\5.2.0\plugins\server
 \eu.echa.iuclid.plugin.resourcesFR.zip

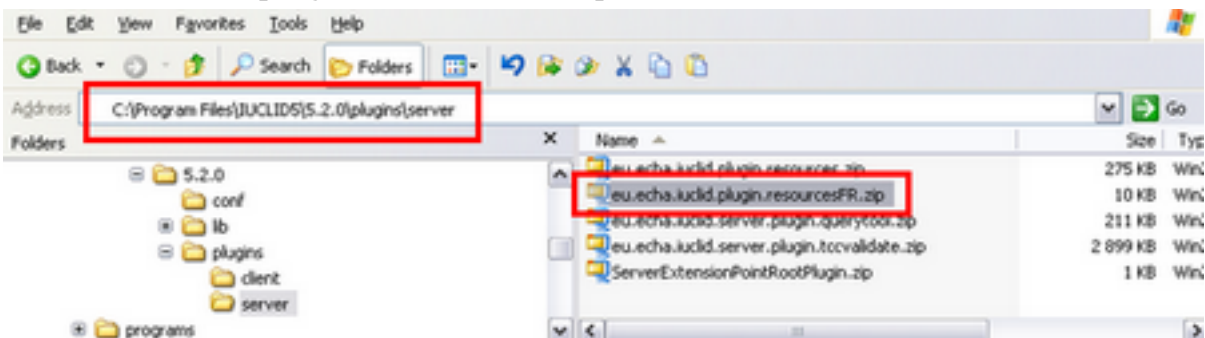

- Copy the file eu.echa.iuclid.plugin.helpsystemFR.zip to the subfolder \plugins\client\ of IUCLID 5.

For example: C:\Program files\IUCLID5\5.2.0\plugins\client
 \eu.echa.iuclid.plugin.helpsystemFR.zip

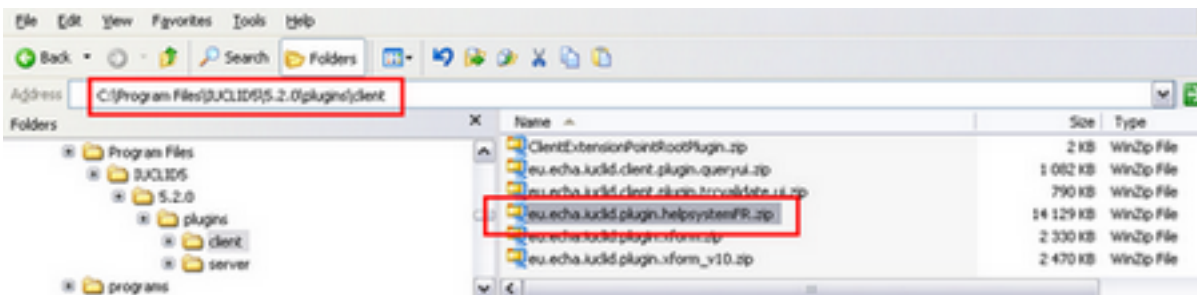



Important

Make sure the zip files are copied as zip files and not extracted.

- Start IUCLID 5. and go to Section 4, “Changing the language setting”.

3.2. Distributed Version


- Save your work and close all IUCLID 5 client application.
- Log in as SuperUser and shut down IUCLID 5 as described in section “20. Preparing Server Shut-Down” of the User manual.
- Locate the installation folder of IUCLID 5 within Tomcat on your hard disk.




Note

The installation folder is usually `C:\Program files\Apache Software Foundation\Tomcat 6.0\webapps\i5server`, however it may be different on your system.

- Copy the file `eu.echa.iuclid.plugin.resourcesFR.zip` to the subfolder `i5server\WEB-INF\classes\plugins` of IUCLID 5.

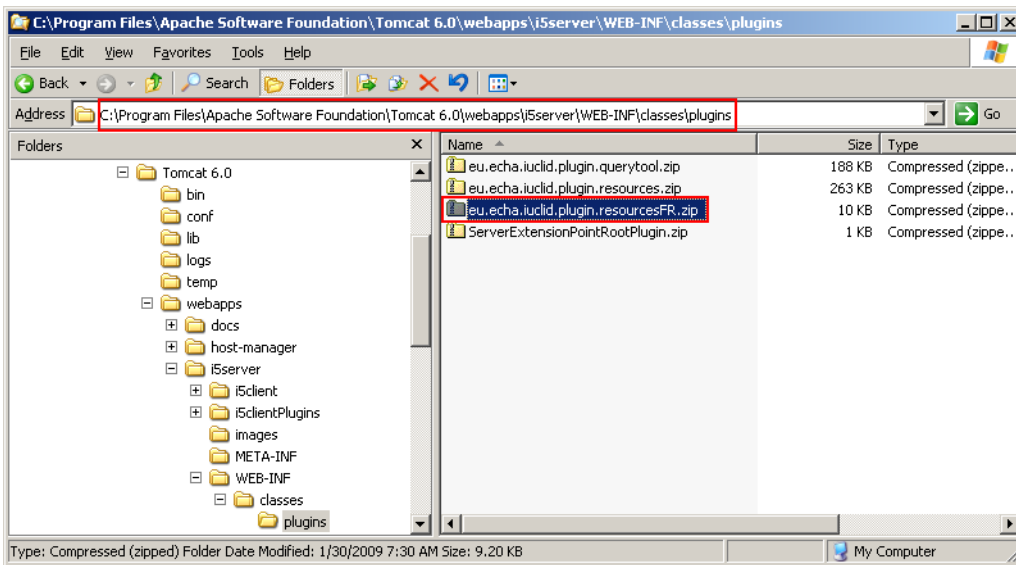

- Copy the file `eu.echa.iuclid.plugin.helpsystemFR.zip` to the subfolder `i5server\i5clientPlugins\remotePlugins` of IUCLID 5.

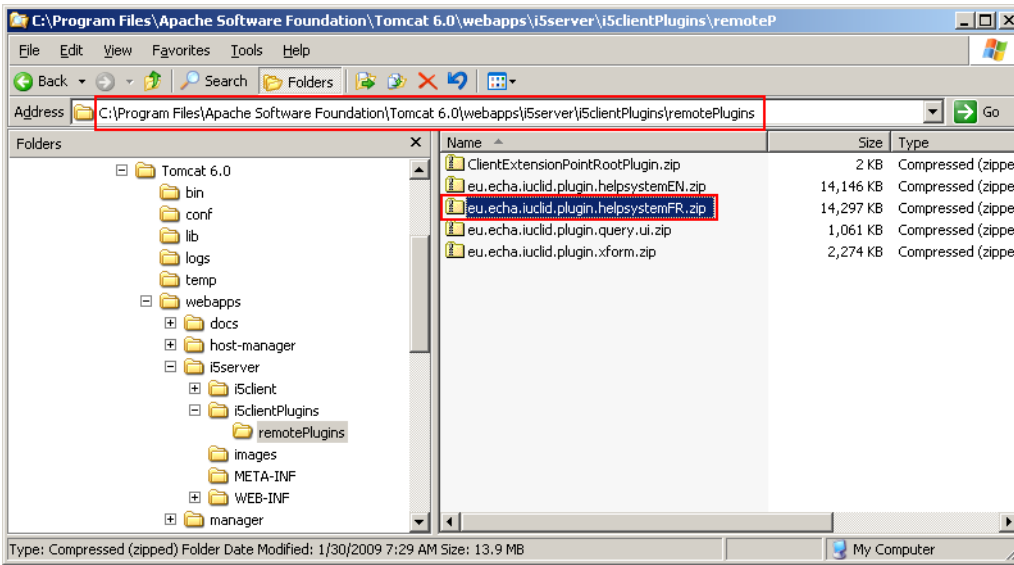


Important

Make sure the zip files are copied as files and not extracted.


- Start IUCLID 5.

4. Changing the language setting


- Start IUCLID 5 and log in to your user account.
- On the main task panel click "user preferences".

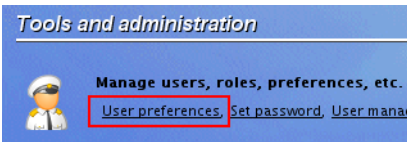

- In the user preferences dialog change the Language setting to French (France).

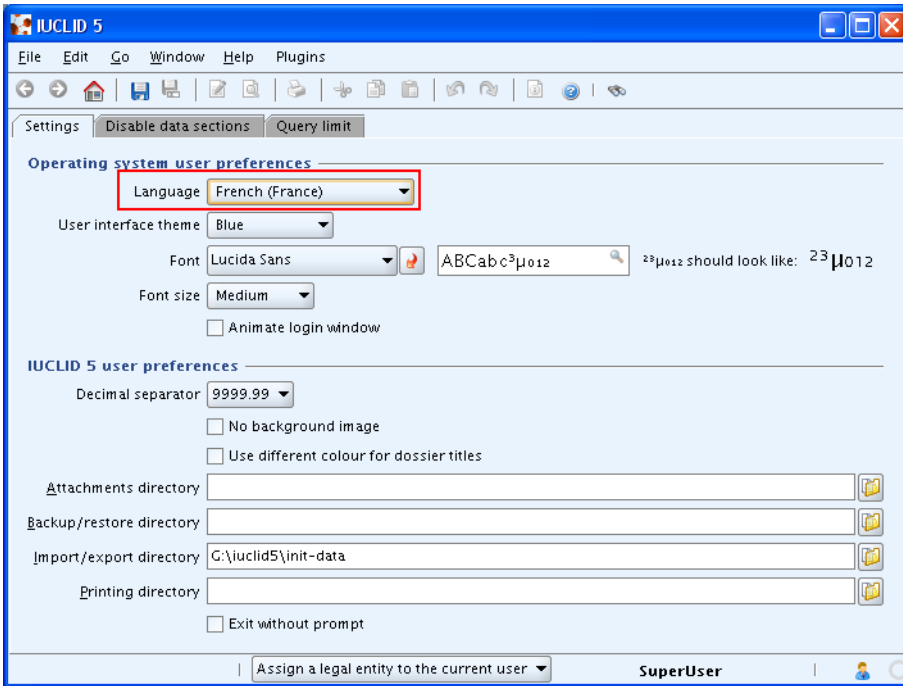

- Close IUCLID 5 and restart in order to activate the new language setting.
- Press the F1 key to visualize the translated helpsystem

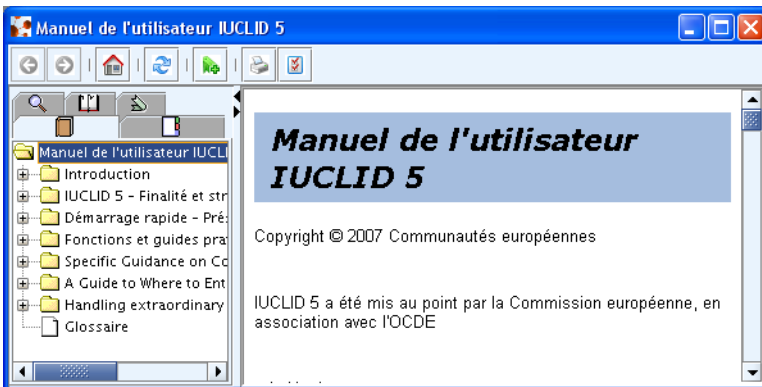

5. Uninstallation

- Change the language setting back to English.
- Save your work and close IUCLID 5.
- Remove the two files installed previously.
- Restart IUCLID 5.

© European Chemicals Agency, 2010. Reproduction is authorised provided the source is fully acknowledged in the form "Source: European Chemicals Agency, <http://echa.europa.eu/>", and provided written notification is given to the ECHA Communication Unit (info@echa.europa.eu).